

Mariner™ - English

Rx Only Caution: Federal Law restricts this device to sale by or on the order of physician

Manufacturer:
SeaSpine Orthopedics Corporation
5770 Armada Drive
Carlsbad, CA 92008, USA

Telephone: +1-760-727-8399

Fax: +1-760-727-8809

Complaints: complaints@seaspine.com

Customer Service: customerservice@seaspine.com

Website: www.seaspine.com


EC REP **European Representative**
mdi Europa GmbH
Langenhagener Str.71
2797 30855 Hannover-Langenhagen, Germany

Telephone: +49 511 39 08 95 30

Fax: +49 511 39 08 95 39

Email: info@mdi-europa.com

Website: www.mdi-europa.com

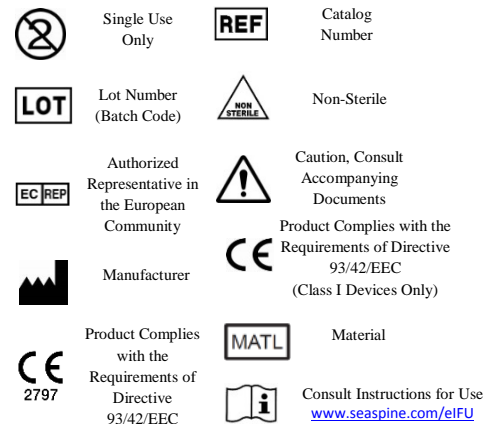

DESCRIPTION

The Mariner Pedicle Screw System is a non-cervical spinal fixation device and instrumentation system intended for use as a posterior pedicle screw fixation system, a posterior non-pedicle screw fixation system, or as an anterolateral fixation system. The system consists of single-use implants including fixed, polyaxial, cephalad/caudal restricted-motion, and medial/lateral restricted motion pedicle screws as well as connecting spinal rods and a separate locking element.

INDICATIONS FOR USE

The intended use of the Mariner Pedicle Screw System is to provide immobilization and stabilization of spinal segments in skeletally mature patients as an adjunct to fusion in the treatment of the following acute and chronic instabilities or deformities of the thoracic, lumbar, and sacral spine.

The indications for use are as follows:

- degenerative disc disease (DDD) as defined by back pain of discogenic origin with degeneration of the disc confirmed by patient history and radiographic studies,
- severe spondylolisthesis (Grades 3 and 4) of the L5-S1 vertebra in skeletally mature patients receiving fusions by autogenous bone graft having implants attached to the lumbar and sacral spine (L3 to sacrum) with removal of the implants after the attainment of a solid fusion,
- spondylolisthesis,
- trauma (i.e., fracture or dislocation),
- spinal stenosis,
- deformities or curvatures (i.e., scoliosis, kyphosis, and/or lordosis),
- spinal tumor,
- pseudarthrosis, and/or
- failed previous fusion.

IMPLANT MATERIALS

Ti-6Al-4V ELI per ASTM F136, Co-28Cr-6Mo per ASTM F1537 and 35Co-35Ni-20Cr-10Mo per ASTM F562.

CONTRAINDICATIONS

Any medical or surgical condition which would preclude the potential benefit of spinal implant surgery is a contraindication. The following conditions may reduce the chance of a successful outcome and should be taken into consideration by the surgeon. This list is not exhaustive:

- **Absolute contraindications:**
 - Infection in or around the operative site
 - Allergy or sensitivity to implant materials
 - Any case not described in the indication
- **Relative contraindications:**
 - Local inflammation
 - Morbid obesity
 - Pregnancy

- Fever or leukocytosis
- Prior fusion at the level(s) to be treated
- Grossly distorted anatomy due to congenital abnormalities
- Rapid joint disease, bone absorption, osteopenia, and/or osteoporosis
- Elevation of sedimentation rate unexplained by other diseases, elevation of white blood count (WBC), or a marked left shift in the WBC differential count
- Any case not requiring bone graft and fusion or where fracture healing is not required
- Patients having inadequate tissue coverage over the operative site or where there is inadequate bone stock, bone quality, or anatomical definition
- Unsuitable or insufficient bone support
- Bone immaturity
- The patient's activity level, mental condition, occupation and/or a patient unwilling to cooperate with the postoperative instructions
- Any case where implant utilization would interfere with anatomical structures or expected physiological performance
- Use of incompatible materials from other systems

POSSIBLE ADVERSE EVENTS

Like other spinal system implants, the following adverse events are possible. This list is not exhaustive:

- Delayed union or nonunion (pseudarthrosis)
- Bending, disassembly or fracture of implant and components
- Loosening of spinal fixation implants may occur due to inadequate initial fixation, latent infection, and/or premature loading, possibly resulting in bone erosion, migration or pain
- Pain, discomfort, or abnormal sensations due to the presence of the device
- Pressure on skin where inadequate tissue coverage exists over the implant, with potential extrusion through the skin.
- Dural leak requiring surgical repair.
- Cessation of growth of the fused portion of the spine.
- Subsidence of the implant into adjacent bone.
- Loss of proper spinal curvature, correction, height and/or reduction.
- Increased biomechanical stress on adjacent levels.
- Improper surgical placement of the implant causing stress shielding of the graft or fusion mass.
- Intraoperative fissure, fracture, or perforation of the spine.
- Postoperative fracture due to trauma, defects, or poor bone stock.
- Serious complications associated with any surgery may occur. These include, but are not limited to: wound complications, infection, genitourinary disorders, gastrointestinal disorders, vascular disorders, including thrombus; bronchopulmonary disorders, including emboli; bursitis, hemorrhage, myocardial infarction, paralysis or death.

WARNINGS AND PRECAUTIONS

Patients with previous spinal surgery at the level(s) to be treated may have different clinical outcomes compared to those without previous surgery. The safety and effectiveness of spinal systems have been established only for spinal conditions with significant mechanical instability or deformity requiring fusion with instrumentation. These conditions are significant mechanical instability or deformity of the spine secondary to severe spondylolisthesis, degenerative spondylolisthesis with objective evidence of neurological impairment, fracture, dislocation, scoliosis, kyphosis, spinal tumor and failed previous fusion (pseudarthrosis). The safety and effectiveness of these devices for any other condition is unknown.

- The implantation of this system should be performed only by experienced spinal surgeons with specific training in the use of this device because this is a technically demanding procedure presenting a risk of serious injury to the patient.
- Based on the fatigue testing results, the surgeon should consider the levels of implantation, patient weight, patient activity level, other patient conditions, etc. which may impact the performance of the system
- Ensure all implants, components or instruments are sterilized prior surgery. The use of non-sterile devices may lead to inflammation, infection or disease.
- Implants should never be reused under any circumstances. A used implant should be discarded. While the implant may appear undamaged, it may have small defects or internal stress patterns and if implanted, could fail to perform as intended and pose safety risks to the patient. The risks include, but are not limited to, mechanical failure, breakage, difficulty with implantation, incompatibility with mating components and infection.

MRI SAFETY

The Mariner Pedicle Screw device has not been evaluated for safety and compatibility in the MR environment. It has not been tested for heating, migration, or image artifact in the MR environment. The safety of the Mariner Pedicle Screw device in the MR environment is unknown. Scanning a patient who has this device may result in patient injury.

CLEANING AND DECONTAMINATION

All instruments and implants that have been previously taken into a sterile surgical field must be decontaminated and cleaned using established hospital methods before sterilization and reintroduction into the sterile surgical field. The following recommendations are for the manual cleaning and decontamination of surgical instruments. These recommendations are considered guidelines with the ultimate responsibility for verifying adequate cleaning remaining with the user. Automated cleaning systems

differ between hospitals and therefore must be qualified by the hospital.

Manual Cleaning Procedure

| | |
|---|--|
| 1 | Remove all gross visible soil with a damp gauze pad or wipe. |
| 2 | Prepare an enzymatic cleaning solution (such as Prolystica® 2X Enzymatic) per manufacturer's instructions. |
| 3 | Immerse the instruments in the cleaning solution and actuate all features so the enzymatic cleaner contacts all mated surfaces and soak for 15 minutes. |
| 4 | Transfer the instruments to fresh cleaning solution (such as Prolystica® 2X Enzymatic). Thoroughly scrub all instruments with a soft bristle cleaning brush while immersed in the enzymatic cleaning solution. Be sure that thorough scrubbing also includes any lumens with an appropriate size brush. Actuate device to allow access to hard to reach areas. |
| 5 | Thoroughly rinse all instruments with warm running water and dry with a clean cloth and/or allow to air dry. |

Automated Cleaning Procedure

| | |
|----|--|
| 1 | Remove all gross visible soil with a damp gauze pad or wipe. Special attention will be required to examine products with tight crevices, voids, and lumens. Lumens may require pre-cleaning with dampened soft bristle brushes and tight crevices, voids, lumens should be flushed with a syringe. |
| 2 | Prepare an enzymatic cleaning solution (such as Prolystica® 2X) per manufacturer's instructions. Immerse the instruments in the cleaning solution and actuate all features so the enzymatic cleaner contacts all mated surfaces and soak for 15 minutes. |
| 3 | Transfer items to a washer and run a cycle with the parameters listed in the following steps. |
| 4 | PRE-WASH: Cold tap water [2 minutes]. |
| 5 | ENZYME WASH: Enzyme wash using cleaner (such as Prolystica® 2X Enzymatic) per manufacturer's recommendations, hot tap water [4 minutes]. |
| 6 | DETERGENT WASH: Detergent wash using detergent (such as Prolystica® 2X Alkaline) per manufacturer's recommendations, hot tap water (66°C/150°F) [2 minutes]. |
| 7 | RINSE 1: Rinse, hot tap water [2 minutes]. |
| 8 | RINSE 2: Purified water rinse (66°C/150°F) [15 seconds]. |
| 9 | DRYING: Hot air dry (82°C/180°F) [12 minutes]. |
| 10 | Remove items from the washer and remove any residual moisture with a lint free clean cloth. |

STERILIZATION

The implants, components and instrumentation in the Mariner System are to be sterilized by the hospital prior to surgery.

The implants and instruments can be sterilized using the provided standard open cases or Aesculap closed cases (standard or PrimeLine lid). Small baskets, trays, and other types of accessories, especially with covers or lids, not provided by SeaSpine for a specific system should not be used. Only SeaSpine standard open cases or Aesculap

closed cases (standard or PrimeLine lid) are validated for use with SeaSpine products.

For information regarding closed cases, please refer to appropriate Instructions for Use provided by the closed case manufacturer.

For standard open cases, trays are to be double wrapped using FDA-cleared sterilization wraps (2 wraps) prior to placement in autoclave. The recommended sterilization cycle will provide a Sterility Assurance Level of (SAL 10⁻⁶).

In a properly functioning and calibrated steam sterilizer, effective sterilization may be achieved using the following parameters:

| | |
|---------------------|---------------|
| Method | Steam |
| Cycle | Pre-vacuum |
| Temperature | 270°F (132°C) |
| Exposure Time | 4 minutes |
| Minimum Drying Time | 30 minutes |

PACKAGING

All packages containing implants should be sealed and intact upon receipt. If the package or product is damaged, the product should not be used and should be returned. The product must be handled, stored, and opened in such a way that it is protected from inadvertent damage or contamination. If a loaner or consignment system is used, all sets should be carefully checked for completeness and all components should be carefully checked for damage before use.

SURGICAL TECHNIQUE

This package insert is designed to assist in using the product and is not intended to provide information on surgical technique. Contact a SeaSpine Representative, customerservice@seaspine.com or 760-727-8399 for a Surgical Technique Guide.

IMPLANT SELECTION

Verify that all parts and necessary instruments are present prior to surgery, including sizes larger and smaller than those that are expected for use. The construct should be assembled prior to surgery.

PREOPERATIVE WARNINGS

- Only patients that meet the criteria described in the indications should be selected.
- Patient condition and/or predispositions such as those described in the contraindications should be avoided.
- Care should be used in the handling and storage of the implants. The implants should not be scratched or damaged. Implants and instruments should be protected during storage and from corrosive environments.
- All non-sterile parts should be cleaned and sterilized before use. Additional sterile components should be available in case of unexpected need.
- Devices should be inspected for damage prior to implantation.
- Care should be used during surgical procedures to prevent damage to the device(s) and injury to the patient.

INTRAOPERATIVE WARNINGS

- Consult Surgical Technique Guide for system specific intraoperative warnings, precautions and recommendations.
- Extreme caution should be used around the spinal cord and nerve roots. Damage to the nerves will cause loss of neurological function.
- Breakage, slippage, or misuse of instruments or implant components may cause injury to the patient or operative personnel.
- Bone graft must be placed in the area to be fused and the graft must be in contact with viable bone.

- Implants and components should not be bent, reshaped, contoured or otherwise modified.
- Use great care to ensure that the implant surfaces are not scratched or notched which may reduce the functional strength of the construct.
- If the construct contains screws, prior to soft tissue closure, recheck all screws to ensure they are tightened. Failure to do so may cause loosening of the other components.

POSTOPERATIVE WARNINGS

- Surgeons should advise patients regarding the risks of surgery and the importance of post-operative compliance.
- The patient should be advised to limit and restrict physical activities, especially lifting and twisting motions and any type of sport participation.
- The patient should be advised that implants may bend, break or loosen despite restriction in activity.
- The patient should be advised to avoid mechanical vibrations that may loosen the device.
- The patient should be advised not to smoke or consume alcohol during recovery.

COMPLAINTS

Immediately notify SeaSpine or a SeaSpine representative by phone, fax or email regarding complaints, malfunctions or adverse events associated with this product. When possible, retain the product involved in the complaint and return to SeaSpine as instructed by SeaSpine Customer Service.

PRODUCT INFORMATION DISCLOSURE

SeaSpine has exercised reasonable care in the selection of materials and the manufacture of these products. SeaSpine warrants to the original purchaser only that each new SeaSpine product is free from manufacturing defects in material and workmanship under normal use and service for a period of six (6) months from the date of delivery by SeaSpine to the original purchaser, but in no event beyond the expiration date stated on any product labeling. SEASPINE DISCLAIMS ALL OTHER WARRANTIES, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). Further, this warranty shall not apply to, and SeaSpine shall not be responsible for, any loss arising in connection with the purchase or use of any SeaSpine product that has been repaired by anyone other than an authorized SeaSpine service representative or altered in any way so as, in SeaSpine's judgment, to affect its stability or reliability, or which has been subject to misuse, negligence or accident, or which has been used otherwise than in accordance with these instructions. IN NO EVENT SHALL SEASPINE BE LIABLE FOR ANY SPECIAL, INCIDENTAL CONSEQUENTIAL, OR CONTINGENT LOSS, DAMAGE, OR EXPENSE, DIRECTLY OR INDIRECTLY ARISING FROM ACQUISITION OR USE OF THIS PRODUCT. SeaSpine neither assumes nor authorizes any person to assume for it any other additional liability or responsibility in connection with this product. SeaSpine intends that this device should be used only by physicians having received proper training in the use of the device.